


ETHICON Megadyne™

REF 2330 IN-LINE FILTER

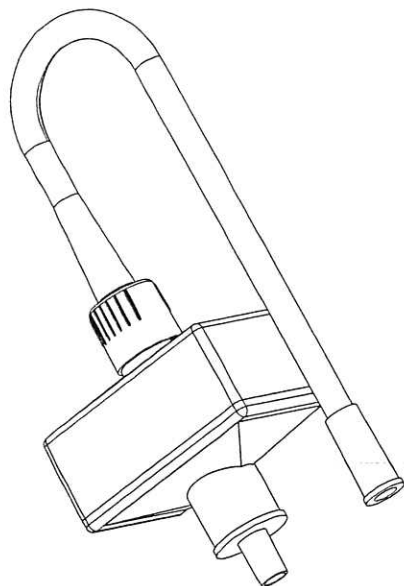

CE REP

Quality First International
Suites 317/318 Barford Business Centre
11 Barford Road
Stratford, London E15 2ST, UK
Tel: +44-208-221-2361
Fax: +44-208-221-1912

3000306-01 REV 002

DRAFT
2017-MKT-17-DCB-094

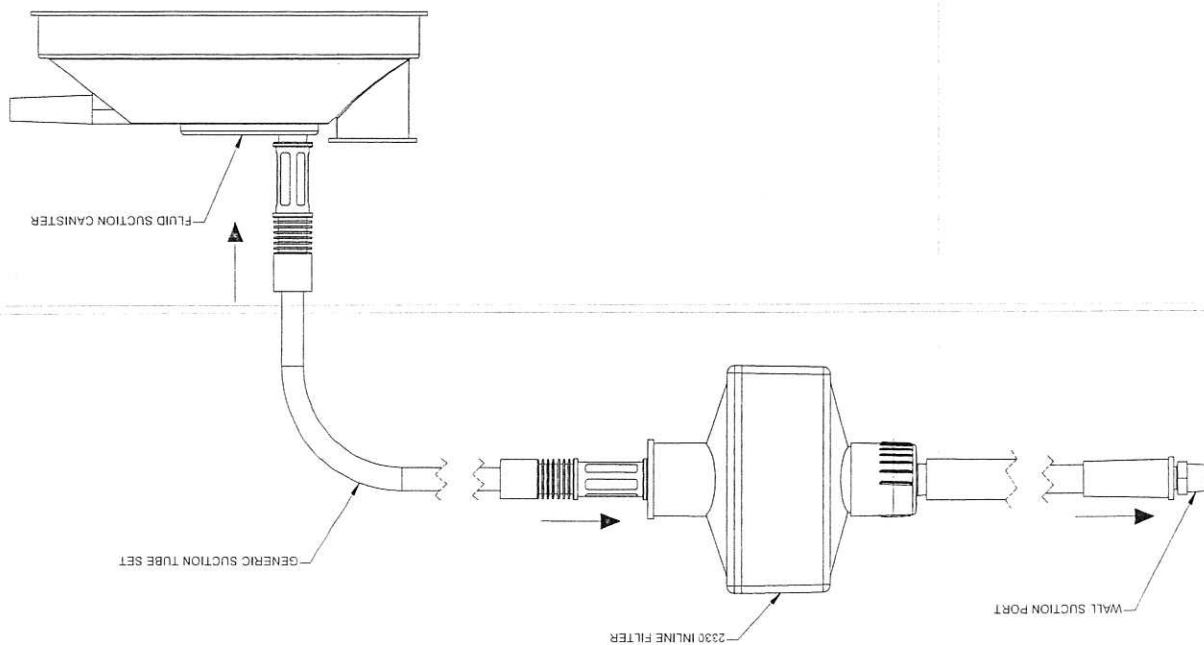


Figure 1

INDICATION FOR USE

The Megadyne In-line Filter is designed to provide ULPV filtration of smoke and aerosols created at a surgical site when used in conjunction with a wall suction system.

INSTRUCTIONS FOR USE:

Insert the In-line Filter between the suction canister and the wall suction connection by:

1. Attach the 2330 In-line Filter's blue connector to the wall suction port.
2. Attach the 2330 In-line Filter and fluid suction canister together via a generic suction tube set. See Figure 1 on back cover.

WARNINGS:




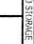





- This filter is considered single-use and should be discarded as per hospital biohazard policy.

CAUTIONS:

- Failure to install a suction canister/fluid collection device may cause the filter to become blocked and reduce filtration.

COMPATIBILITY:

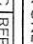- Megadyne recommends the use of the In-line Filter with suction canister when no other smoke removal system is present.

	Do not reuse
	Do not use if package is damaged
	Consult instructions for use
	Normal storage conditions are assumed. Brief excursions to limits allowed
	Quantity of individual units per labeled container
	Batch code
Rx ONLY	Federal (USN) law restricts this device to sale by or on the order of a physician
	Authorized representative in the European community
	Manufacturer
	Catalog number